

The Dräger Vapor 3000 / D-Vapor 3000 Advanced vaporizer with smart functions

The new Vapor 3000 vaporizers offer advanced features to better prepare you for the challenges of today's OR – such as visibility in darker environments. In combination with the Perseus® A500 anesthesia workstation, the Vapor 3000 vaporizers offer a new level of confidence.



NEW FEATURES FOR ADDED COMFORT


The new Vapor 3000 series incorporates features from previous models, such as transport functionality and a large filling volume of 300 ml. However, it now includes advanced features that help you perform anesthesia with more confidence than ever before. In addition to ensuring secure storage and application of relevant current volatile agents, it is now equipped with an easy-open sealing cap, as well as setting and filling level illumination. The control dial features an off (0) setting and a transport setting (T). Additionally, the Vapor 3000 series can now even communicate with the new Perseus A500 anesthesia workstation.

COMMUNICATION WITH THE PERSEUS A500

Optimized to seamlessly communicate with the Perseus A500, the Vapor 3000 series offers a new level of convenience. When mounted on the Perseus A500 with the Vapor Communication function, vaporizers are automatically recognized and identified. The type and control dial setting of each vaporizer can easily be viewed directly on the screen of the Perseus A500. Now, clinicians have a single location to view vaporizer settings, measured anesthetic agent concentrations and other essential ventilation information – without sacrificing the feel and functionality of mechanical vaporizers. Additionally, low vaporizer filling levels are automatically indicated by an acoustic signal and visual alarms on the Perseus A500 screen and vaporizers.



Communication with the Perseus® A500
Vaporizer identification and the control dial position of 2.4.

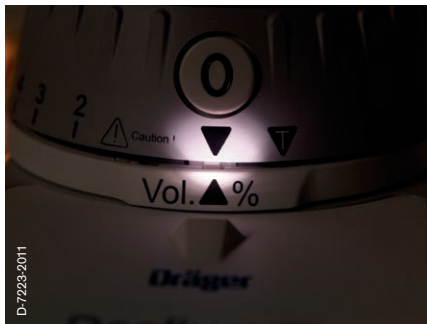

ILLUMINATED CONCENTRATION SETTING AND FILLING LEVEL INDICATOR

The Vapor 3000 series comes with integrated illumination for both the control dial and filling level indicators. This means you can be sure of your supply and your dosage levels even in dark environments and without the need for additional light sources. The level indicator now magnifies the agent level for better, wide-angle visibility. A blue stripe at the back of the glass further enhances the view of the agent level.

NO RECALIBRATION

The new family of vaporizers helps you save time: they do not need to be recalibrated at any time during their period of use.

Improved sealing-cap design and operation

The newly designed sealing cap can be opened and sealed with just a 90° turn, while its grip clearly indicates the status of the seal; the horizontal grip position tells you that the cap is open, the vertical position means that the cap is sealed. This means you can be sure of the status of the seal at a glance.



CORPORATE HEADQUARTERS

Drägerwerk AG & Co. KGaA
Moislinger Allee 53–55
23558 Lübeck, Germany

www.draeger.com

Manufacturer:

Dräger Medical GmbH
Moislinger Allee 53–55
23558 Lübeck, Germany
The quality management system at Dräger Medical GmbH is certified according to ISO 13485, ISO 9001 and Annex II.3 of Directive 93/42/EEC (Medical devices).

USA

Draeger Medical, Inc.
3135 Quarry Road
Telford, PA 18969-1042, USA
Tel +1 215 721 5400
Toll-free +1 800 437 2437
Fax +1 215 723 5935
info.usa@draeger.com